




The Biotechnology Education Company®

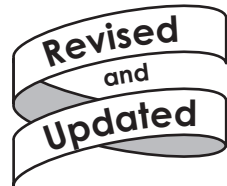

EDVO-Kit #
331

Cloning of a PCR Amplified Gene

**Storage: See Page 3 for
specific storage instructions**

EXPERIMENT OBJECTIVE:

The objective of this experiment is to amplify a DNA fragment by Polymerase Chain Reaction (PCR) and to clone the amplified DNA by using the blue/white cloning system.

This experiment is designed for DNA staining with InstaStain® Ethidium Bromide.

EDVOTEK, Inc. • 1-800-EDVOTEK • www.edvotek.com

Table of Contents

Experiment Components	3
Experiment Requirements	4
Background Information	5
Experiment Procedures	
Experiment Overview and General Instructions	12
Laboratory Safety	13
Module I: Amplification of a 270 bp DNA Fragment by PCR and Separation by Electrophoresis	14
Module II: Preparation of Insert for Ligation	17
Module III: Ligation of the PCR Amplified 270 bp DNA into pUC19	19
Module IV: Transformation	20
Study Questions	22
Instructor's Guidelines	
Notes to the Instructor	24
Pre-Lab Preparations	27
Experiment Results and Analysis	32
Study Questions and Answers	33
Appendices	
A PCR Experimental Success Guidelines	36
B Polymerase Chain Reaction Using Three Waterbaths	37
C Preparation and Handling of PCR Samples With Wax	38
D 1.0% Agarose Gel Preparation	39
E 1.0% Agarose Gels - Quantity Preparations	40
F Staining and Visualization of DNA with InstaStain® Ethidium Bromide Cards	41
• Material Safety Data Sheets	42



EVT 003157K

EDVOTEK, The Biotechnology Education Company, and InstaStain are registered trademarks of EDVOTEK, Inc.. Ready-to-Load and UltraSpec-Agarose are trademarks of EDVOTEK, Inc.

Components & Requirements

This experiment contains reagents to perform five PCR reactions and five cloning experiments.

Sample volumes are very small. For liquid samples, it is important to quickly spin the tube contents in a microcentrifuge to obtain sufficient volume for pipeting. Spin samples for 10-20 seconds at maximum speed.

All components are intended for educational research only. They are not to be used for diagnostic or drug purposes, nor administered to or consumed by humans or animals.

THIS EXPERIMENT DOES NOT CONTAIN HUMAN DNA. None of the experiment components are derived from human sources.

Reagents for PCR

		Storage
P1	DNA Template for Amplification	-20°C
P2	Primer Set (two primers)	-20°C
P3	Tubes with PCR Reaction Pellets™ (Each PCR reaction pellet™ contains dNTP Mixture, <i>Taq</i> DNA Polymerase Buffer, <i>Taq</i> DNA Polymerase and MgCl ₂)	Room temp.
P4	Enzyme Grade UltraPure Water	-20°C
P5	200 bp DNA ladder	-20°C
	• Wax Beads (for thermal cyclers without a heated lid)	Room temp.
	• 10x Gel Loading Solution	Room temp.

Reagents to Clean and Prepare DNA for Ligation

C1	Purification Buffer	Room temp.
C2	Wash Buffer	Room temp.
C3	Restriction Enzyme Reaction Buffer	-20°C
C4	Enzyme Grade UltraPure Water	-20°C
C5	<i>Eco</i> RI	-20°C
C6	<i>Hind</i> III	-20°C
C7	Restriction Enzyme Dilution Buffer	-20°C
	• Spin Columns and Reservoirs	Room temp.

Reagents for Ligation

L1	pUC19 Plasmid Vector cut with <i>Eco</i> RI and <i>Hind</i> III	-20°C
L2	DNA Ligase	-20°C
L3	Ligation buffer	-20°C

Reagents for Transformation

T1	Ampicillin	-20°C
T2	IPTG	-20°C
T3	X-Gal in Solvent	-20°C
T4	CaCl ₂	-20°C
1	Bottle of Recovery Broth	Room temp.
1	Vial of Bacterial LyphoCells™	Room temp.
1	Bottle of Ready Pour Agar	Room temp.

Other Components:

- UltraSpec-Agarose™
- Electrophoresis Buffer (50x)
- InstaStain® Ethidium Bromide
- Microcentrifuge Tubes
- PCR tubes (0.2 ml - for thermal cyclers with 0.2 ml template)
- Petri plates

Cloning of a PCR Amplified Gene


Requirements

- Thermal cycler (EDVOTEK Cat. # 532 highly recommended) or three waterbaths
- Horizontal gel electrophoresis apparatus
- D.C. power supply
- Balance
- Microcentrifuge
- Water bath (37°C & 65°C)
- Incubation oven (37°C)
- UV Transilluminator or UV Photodocumentation system
- UV safety goggles
- Automatic micropipets (5-50 μ l & 0.5-10 μ l) with tips
- Microwave, hot plate or burner
- Pipet pump
- 250 ml flasks or beakers
- Hot gloves
- Disposable vinyl or latex laboratory gloves
- Ice buckets and ice
- Distilled or deionized water
- Isopropanol

Online Ordering
now available

www.edvotek.com

Visit our web site for information about EDVOTEK's complete line of "hands-on" experiments for biotechnology and biology education.

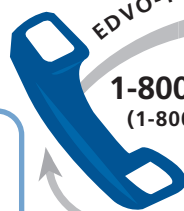

EDVO-TECH SERVICE

1-800-EDVOTEK
(1-800-338-6835)

Mon - Fri 9 am - 6 pm ET

**Technical Service
Department**

Mon - Fri
9:00 am to 6:00 pm ET

FAX: (301) 340-0582
Web: www.edvotek.com
email: edvotek@aol.com

**Please have the following
information ready:**

- Experiment number and title
- Kit lot number on box or tube
- Literature version number (in lower right corner)
- Approximate purchase date



EVT 003157K

Polymerase Chain Reaction, Ligation and Cloning

The Polymerase Chain Reaction (PCR) has made an extraordinary impact on various aspects of biotechnology. The success of utilizing PCR is due to the specificity endowed in the hybridization of nucleic acid and simplicity of the reaction. PCR has also made amplification of genes an alternate approach to traditional cloning experiments. It is currently being used in genome projects, in DNA mapping and sequencing. PCR is also applied in DNA based forensics, paternity and determination of evolutionary relationships.

In a typical PCR reaction, the first step is the preparation of the DNA sample that is extracted from various biological sources. Segments of DNA or genes to be amplified are referred to as the targets and the two synthetic oligonucleotides required for the PCR reaction are referred to as primers. Because PCR is very sensitive, only a few copies of the gene is required. Freshly isolated DNA will give the best amplification results compared to DNA extracted from older specimens that may be degraded. The set of two primers, usually in the range between 15 and 30 nucleotides, are chemically synthesized to correspond to the two ends of the gene or DNA to be amplified. The primer concentrations are always in excess of the DNA target. The nucleotide primer sequences for the DNA amplification reaction are determined to yield the best hybridization.

ABOUT DNA POLYMERASES

DNA Polymerase I is the best studied DNA polymerase. With some exceptions, the mechanism of DNA polymerization is essentially identical for all polymerases. This DNA polymerase is a single polypeptide with a molecular weight of 109,000 (approximately 1000 amino acids). Like other DNA polymerases it requires a primed template, the four deoxynucleotide triphosphates and magnesium for DNA synthesis. Polymerization of the nucleotides occurs in the 5' to 3' direction by the addition of a 5' phosphorylated nucleotide to the free 3' hydroxyl group of the growing DNA chain with the concomitant production of pyrophosphate. The primer is antiparallel and is base paired to the template strand. In vivo, RNA fragments containing about 10 nucleotides serve to prime DNA synthesis. RNA primers are synthesized by the enzyme primase which is a specialized DNA dependent RNA polymerase.

DNA polymerase I also has a 3'-5' exonucleolytic activity that cleaves DNA starting with a free terminal 3' hydroxyl group on the growing chain to yields 5'-deoxynucleotide monophosphates. This activity serves as a proof-reading function during DNA synthesis. It recognizes distortions in the 3' ends of the growing DNA chain caused by mismatched bases between the template and the growing chain. The enzyme also possesses a 5' to 3' nuclease activity. This activity is part of a repair enzyme system that removes damaged DNA bases.

Polymerase Chain Reaction, Ligation and Cloning

Limited proteolysis of DNA polymerase I by subtilisin or trypsin produces two polypeptide fragments having molecular weights of 76,000 and 36,000. The larger polypeptide, known as the Klenow fragment contains the polymerization and 3' to 5' exonuclease activities while the smaller fragment contains the 5' to 3' nuclease activity. The Klenow fragment was used in the initial PCR experiments and was subsequently replaced by *Taq* DNA polymerase that is thermally stable. The source of *Taq* DNA polymerase is *Thermus aquaticus*. Several other thermo-stable polymerases have also been purified from other thermophilic bacteria. Because of its thermal stability, *Taq* DNA polymerase will not be denatured during the high temperature PCR steps described below. In the initial experiments, the Klenow fragment had to be replenished after each DNA denaturation step (94°C). The thermostability of *Taq* DNA polymerase eliminates replenishment, thereby simplifying PCR.

ABOUT THE POLYMERASE CHAIN REACTION

A typical PCR reaction mixture contains DNA, the four deoxynucleotide triphosphates, Mg+2 *Taq* DNA polymerase and the reaction buffer. The total incubation reaction is usually small (10 to 25 µl) in volume. If water baths or a thermal cycler without a heated lid is used, the incubation reaction mixture is layered with mineral oil to minimize evaporation. With thermal cyclers that have heated lids layering with mineral oil is not required.

The PCR reaction mixture is then exposed to sets of three step temperature cycles. The first temperature 94°C melts the hydrogen bonds between the two DNA strands. The temperature is then reduced between 45°C to 60°C to hybridize the two primers to each of target DNA strands. The temperature is then increased to 72°C, the optimum for *Taq* DNA polymerase. During this step DNA is synthesized. The three temperature steps of a cycle are usually repeated 20 to 30 times (Figure1). This process is made efficient by placing the reaction tubes in thermal cyclers that are programmed to alternate and maintain temperatures. Amplified DNA products are detected by gel electrophoresis analysis.

It should be noted that PCR amplification can introduce a small number of mutations in the amplified DNA product. Mutations introduced can be controlled by the use of low nucleotide concentrations. Various undesired amplification products such as short DNA fragments are also synthesized. For best results in subsequent steps such as ligation in plasmids, the amplified DNA is purified from unused dNTPs, primers and *Taq* DNA polymerase prior to further use.

The correct PCR product will be in three forms. These products are due to an intrinsic property of the *Taq* DNA polymerase that appends an extra nucleotide usually a dA, at the 3' ends of blunt double-stranded DNA. This will yield three forms where the first form will be double-stranded with no dA appended at either 3' end. The second form will have a dA on one of the



Duplication of this document, in conjunction with use of accompanying reagents, is permitted for classroom/laboratory use only. This document, or any part, may not be reproduced or distributed for any other purpose without the written consent of EDVOTEK, Inc. Copyright © 1996, 1997, 1998, 2000, 2001, 2002, 2005, 2007 EDVOTEK, Inc., all rights reserved
EVT 003157K

Polymerase Chain Reaction, Ligation and Cloning

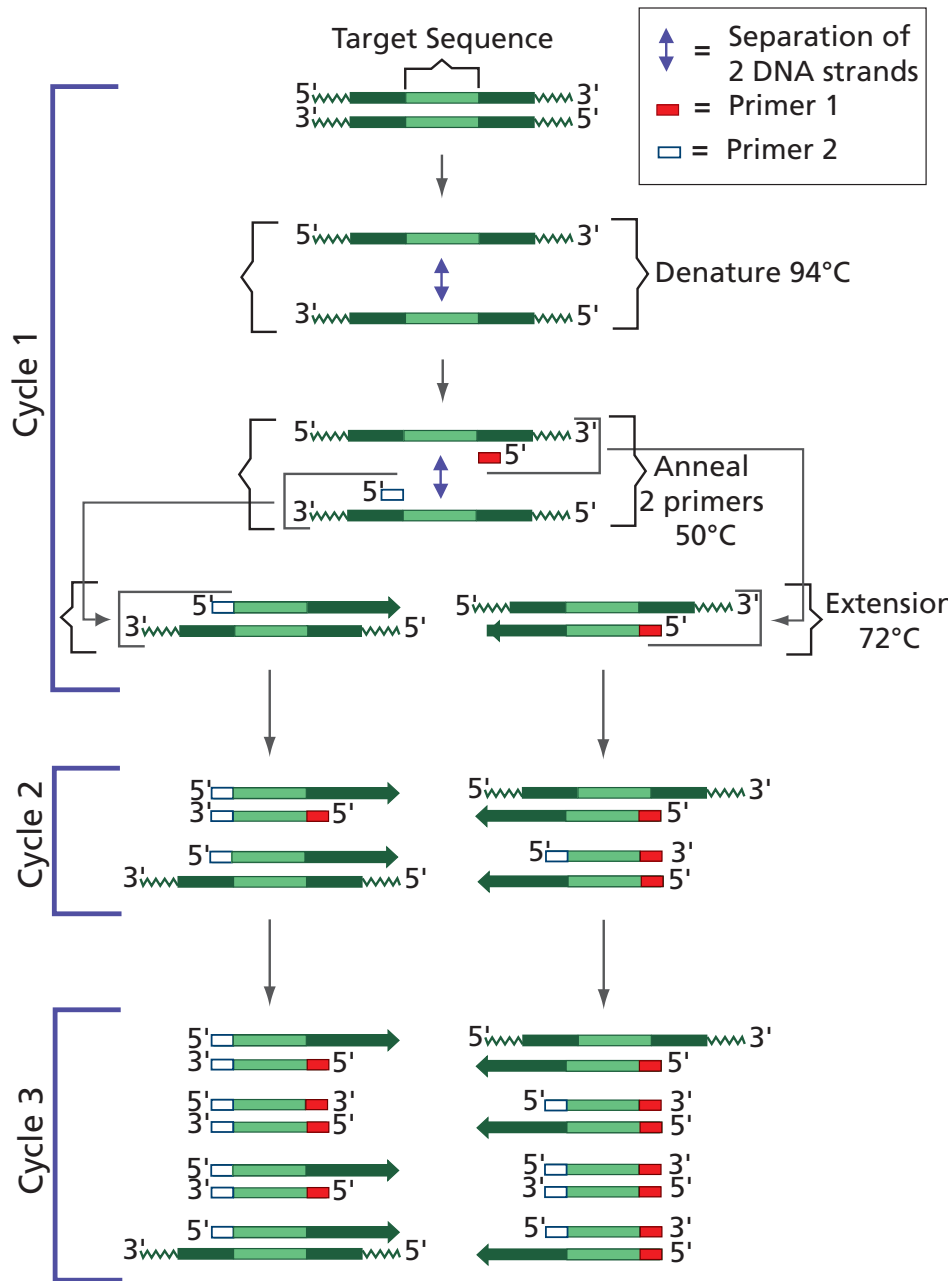


Figure I: The Polymerase Chain Reaction

Background Information

Duplication of this document, in conjunction with use of accompanying reagents, is permitted for classroom/laboratory use only. This document, or any part, may not be reproduced or distributed for any other purpose without the written consent of EDVOTEK, Inc. Copyright © 1996, 1997, 1998, 2000, 2001, 2002, 2005, 2007 EDVOTEK, Inc., all rights reserved



Polymerase Chain Reaction, Ligation and Cloning

two DNA strands and the third form will have a dA residues on both DNA strands. This single nucleotide addition is referred to as “template independent polymerization”. The presence of this nucleotide will have a negative impact on a subsequent cloning step. When overhanging end cloning (sticky ends) is desired, the best strategy would be to use primers that flank restriction enzyme sites present at the two ends of the target DNA that can be amplified. The amplified product is then digested by a restriction enzyme.

ABOUT PLASMIDS

The plasmid pUC19 used for this experiment is derived from the pUC series. It has a single recognition site for *Eco* RI and *Hind* III (restriction enzymes), that are located in a polylinker. The polylinker region known as MCR (multiple cloning region) contains several

single restriction enzyme sites that facilitates the insertion of DNA. The pUC19 plasmid (Figure 2) is present in multiple copies in a host *E. coli* cell, and has been cleverly modified by genetic engineering. These modifications include the addition of the lac Z gene that codes for beta-galactosidase, an enzyme involved in lactose metabolism. DNA inserted into the MCR interrupts the lac Z' gene and prevents the formation of a functional beta-galactosidase protein. As a result clones of interest will appear as white colonies instead of blue on selection agar plates that contain ampicillin.

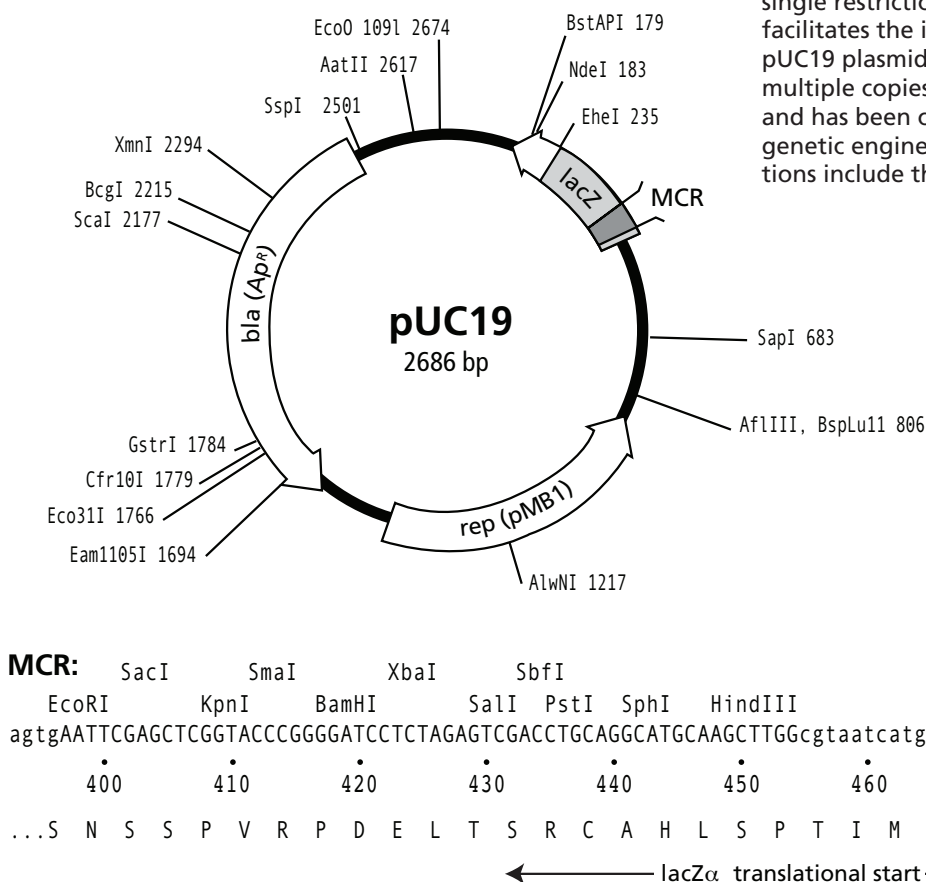


Figure 2: Plasmid pUC19 and MCR Map



Duplication of this document, in conjunction with use of accompanying reagents, is permitted for classroom/laboratory use only. This document, or any part, may not be reproduced or distributed for any other purpose without the written consent of EDVOTEK, Inc. Copyright © 1996, 1997, 1998, 2000, 2001, 2002, 2005, 2007 EDVOTEK, Inc., all rights reserved
EVT 003157K

Polymerase Chain Reaction, Ligation and Cloning

CONSTRUCTION OF A RECOMBINANT PLASMID

Ligation of the PCR amplified DNA to the linearized plasmid is catalyzed by T4 DNA ligase. The enzyme catalyzes the formation of phosphodiester bonds by the condensation of a 5' phosphate and 3' hydroxyl group of adjacent nucleotides. Each phosphodiester bond formation results in the hydrolysis of ATP to AMP and pyrophosphate. Ligation of DNA fragments having cohesive termini is usually achieved at temperatures between 4°C to 22°C. These temperatures will allow for annealing between complementary DNA ends that serves as a prerequisite for ligation.

When the plasmid and insert have the same cohesive termini, the orientation of the sub-cloned DNA (ligated in the plasmid) will vary due to the symmetrical nature of the termini. Statistically one would expect to find a 50:50 occurrence for the DNA orientation in bacterial colonies obtained from the same transformation reaction. Therefore the insert in the recombinant plasmid can be in either one of two directions relative to a fixed point in the vector (Figure 3A).

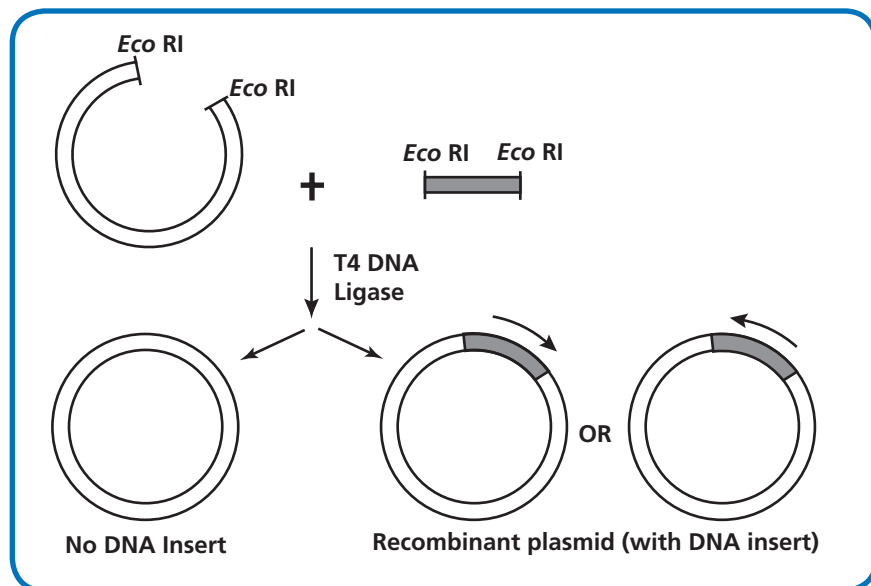


Figure 3A: Use of a single restriction enzyme (*Eco* RI) for the ligation of DNA in an *Eco* RI pre-digested plasmid. As noted by the two arrows the DNA insert is ligated bi-directionally.

Polymerase Chain Reaction, Ligation and Cloning

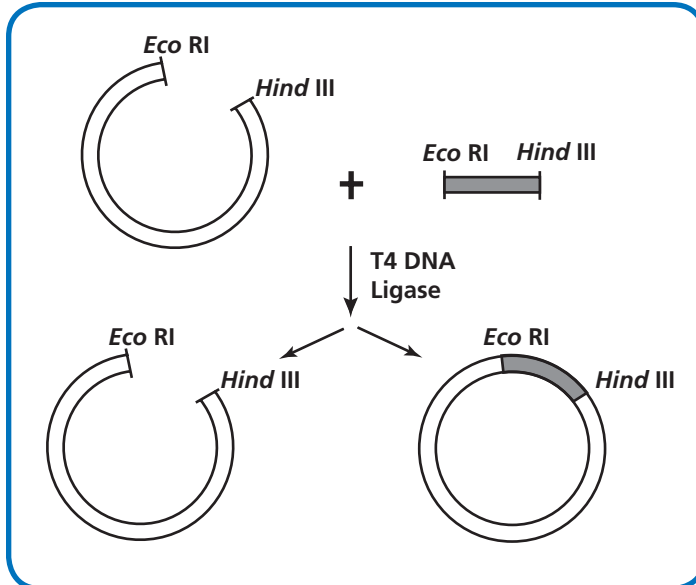


Figure 3B: Demonstrates *Eco RI* and *Hind III* co-digested plasmid and DNA fragment that upon ligation will yield the desired recombinant.

In the simplest reaction, ligation of a plasmid and PCR amplified DNA (cleaved by the same restriction enzyme) will form a circular recombinant plasmid (Figure 3A). The stoichiometry of this reaction is complex and is based on the length and relative concentrations of the two DNAs, the amount enzyme and the ionic strength of the reaction. In this reaction the plasmid (without the amplified DNA) will also circularize (Figure 3A).

To obtain a recombinant with a DNA insert in the desired direction the plasmid is co-digested within the multiple cloning region by two restriction enzymes that will produce cohesive ends that do not match. A similarly co-digested DNA fragment will be ligated in the desired orientation (Figure 3B). In this reaction the plasmid will not circularize because the ends are derived from the co-digestion by two different restriction enzymes. White colonies that may be present would be due to undigested plasmid.

ABOUT TRANSFORMATION:

For this experiment competent cells are prepared from cultures of *E. coli*, strain JM109. This strain does not have any natural antibiotic resistance or plasmids and lacks restriction enzymes. Transformation with the recombinant DNA allows for its expression, propagation, and purification. Linear plasmids and large concatamers do not transform competent cells, while supercoiled DNA has the highest transformation efficiencies. Only small amounts of DNA, typically less than 10 nanograms, are required for transformation. In fact, transformation is inhibited by DNA exceeding 100 nanograms. Even with this amount of DNA, only 1 in 10,000 cells successfully incorporate the recombinant DNA.

Transformation efficiency is based on the number of transformants obtained per microgram of DNA. As an example to determine transformation efficiency, 10 nanograms of DNA were used for a reaction and cells were allowed to recover in a final volume of 1 ml but only one tenth of this volume



Polymerase Chain Reaction, Ligation and Cloning

was plated and produced 100 colonies on a selective agar medium. Therefore, 1000 transformants are present per ml. Keeping in mind that each colony grew from one transformed cell, the efficiency would be $1000/0.01\mu\text{g} = 1 \times 10^5$. Transformation efficiencies of 10^5 to 10^6 are sufficient for most classroom cloning experiments. When cloning of single copy genes from genomic DNA is done, required efficiencies are 10^7 to 10^8 .

This experiment has four modules with the following objectives:

1. To amplify a DNA fragment using the PCR reaction.
2. To prepare PCR amplified DNA for ligation.
3. To ligate the PCR amplified DNA in pUC19
4. To transform *E. coli* host cells and determine the number of white colonies (recombinant plasmid).

Experiment Overview and General Instructions

BEFORE YOU START THE EXPERIMENT

1. Read all instructions before starting the experiment.
2. If you will be conducting PCR using a thermal cycler without a heated lid, also read the Appendix entitled "Preparation and Handling PCR Samples with Wax".
3. If you will be using three waterbaths to conduct PCR, read the two appendices entitled "Polymerase Chain Reaction Using Three Waterbaths" and "Handling samples with wax overlays".
4. Write a hypothesis that reflects the experiment and predict experimental outcomes.

EXPERIMENT OBJECTIVE:

The objective of this experiment is to amplify a DNA fragment by Polymerase Chain Reaction (PCR) and to clone the amplified DNA by using the blue/white cloning system.

BRIEF DESCRIPTION OF EXPERIMENT:

This experiment has four modules with the following objectives:

1. To amplify a DNA fragment using the PCR reaction.
2. To prepare PCR amplified DNA for ligation.
3. To ligate the PCR amplified DNA in pUC19.
4. To transform *E. coli* host cells and to quantitate the number of white colonies.




Duplication of this document, in conjunction with use of accompanying reagents, is permitted for classroom/laboratory use only. This document, or any part, may not be reproduced or distributed for any other purpose without the written consent of EDVOTEK, Inc. Copyright © 1996, 1997, 1998, 2000, 2001, 2002, 2005, 2007 EDVOTEK, Inc., all rights reserved


EVT 003157K

Laboratory Safety

1. Gloves and goggles should be worn routinely as good laboratory practice.
2. Exercise extreme caution when working with equipment that is used in conjunction with the heating and/or melting of reagents.
3. **DO NOT MOUTH PIPET REAGENTS - USE PIPET PUMPS.**
4. Exercise caution when using any electrical equipment in the laboratory.
 - Although electrical current from the power source is automatically disrupted when the cover is removed from the apparatus, first turn off the power, then unplug the power source before disconnecting the leads and removing the cover.
 - Turn off power and unplug the equipment when not in use.
5. EDVOTEK injection-molded electrophoresis units do not have glued junctions that can develop potential leaks. However, in the unlikely event that a leak develops in any electrophoresis apparatus you are using, **IMMEDIATELY SHUT OFF POWER.** Do not use the apparatus.
6. Always wash hands thoroughly with soap and water after handling reagents or biological materials in the laboratory.

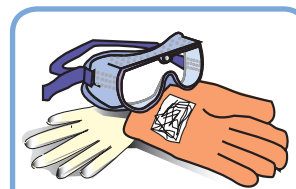

Wear gloves
and safety goggles

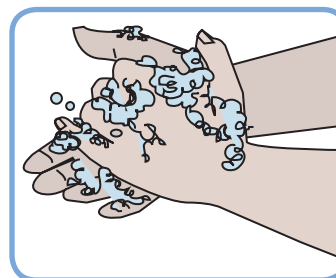

Module I: Amplification by PCR and Separation by Electrophoresis

The PCR reaction pellet™ contains Taq DNA polymerase, the four deoxytriphosphates, Mg⁺² and buffer.

Sample volumes are very small. For liquid samples, it is important to quickly spin the tube contents in a microcentrifuge to obtain sufficient volume for pipeting. Spin samples for 10-20 seconds at maximum speed.

AMPLIFICATION OF A 270 BP DNA FRAGMENT BY PCR

1. Program thermal cycler for 30 cycles with the following schedule:

<u>Initial Denaturation</u>	<u>30 cycles @</u>	<u>Final Extension</u>
94°C for 5 min.	94°C for 30 sec.	72°C for 5 min.
	50°C for 30 sec.	
	72°C for 30 sec.	

2. Obtain one of the tubes (P3) containing the PCR Reaction pellet™ and label the side of the tube and top of the cap with your lab group number or initials.
3. Add the following to tube P3:

5 µl	DNA Template for Amplification (P1)
10 µl	Primer Set (P2)
15 µl	Enzyme Grade Ultrapure Water (P4).
4. Gently mix the reaction tube and pulse spin it in a microcentrifuge to collect the entire sample at the bottom of the tube.
5. If your thermal cycler has a 0.5 ml template, proceed to step 6.
If your thermal cycler has a 0.2 ml template, transfer the entire contents of your PCR reaction tube into a clean 0.2 ml PCR tube before proceeding to step 6.
6. Check to see if your thermal cycler is equipped with a heated lid. If it is, proceed directly to placing the tube in the thermal cycler.
If your thermal cycler **does not** have a heated lid, add one wax bead to the tube before placing the tube in the thermal cycler.
7. After the tubes are placed in the thermal cycler, process the samples for 30 cycles according to the schedule outlined in step 1.
8. To a clean tube add:

3 µl	10x Gel Loading Solution
12 µl	Distilled water.
9. After the final PCR cycle, remove 5 µl of the amplified DNA and transfer it to the tube containing the 10x gel loading solution and water (from step 8). This sample can be analyzed on a 1.0% agarose gel. The remaining amplified DNA (insert) is ready to be cleaned and purified for ligation in Module II.



OPTIONAL STOPPING POINT

Freeze at -20°C the amplified DNA and sample to be analyzed on a 1% gel, or continue with electrophoresis and/or Preparation of Insert for Ligation (Module II).



Duplication of this document, in conjunction with use of accompanying reagents, is permitted for classroom/laboratory use only. This document, or any part, may not be reproduced or distributed for any other purpose without the written consent of EDVOTEK, Inc. Copyright © 1996, 1997, 1998, 2000, 2001, 2002, 2005, 2007 EDVOTEK, Inc., all rights reserved
EVT 003157K

Module I: Amplification by PCR and Separation by Electrophoresis

If you are unfamiliar with agarose gel preparation and electrophoresis, detailed instructions and helpful resources are available at www.edvotek.com

AGAROSE GEL REQUIREMENTS

- Recommended gel size: 7 x 7 cm
- Placement of well-former template: first set of notches
- Agarose gel concentration: 1.0%

PREPARING THE AGAROSE GEL

1. Close off the open ends of a clean and dry gel bed (casting tray) by using rubber dams or tape.
2. Place a well-former template (comb) in the first set of notches at the end of the bed. Make sure the comb sits firmly and evenly across the bed.
3. To a 250 ml flask or beaker, add agarose powder and buffer as indicated in the Reference Tables (Appendix A) provided by your instructor. Swirl the mixture to disperse clumps of agarose powder.
4. With a marking pen, indicate the level of the solution volume on the outside of the flask.
5. Heat the mixture using a microwave oven or burner to dissolve the agarose powder.
6. Cool the agarose solution to 60°C with careful swirling to promote even dissipation of heat. If detectable evaporation has occurred, add distilled water to bring the solution up to the original volume marked in step 4.

Important Note



Continue heating until the final solution appears clear (like water) without any undissolved particles. Check the solution carefully. If you see "crystal" particles, the agarose is not completely dissolved.

After the gel is cooled to 60°C:


7. Place the bed on a level surface and pour the cooled agarose solution into the bed.
8. Allow the gel to completely solidify. It will become firm and cool to the touch after approximately 20 minutes.
9. After the gel is solidified, be careful not to damage or tear the wells while removing the rubber dams or tape and comb(s) from the gel bed.
10. Place the gel (on its bed) into the electrophoresis chamber, properly oriented, centered and level on the platform.
11. Fill the electrophoresis apparatus chamber with the appropriate amount of diluted (1x) electrophoresis buffer (refer to Table B on the instruction Appendix provided by your instructor).

Module I: Amplification by PCR and Separation by Electrophoresis

The Experiment

Reminder:

Before loading the samples, make sure the gel is properly oriented in the apparatus chamber.

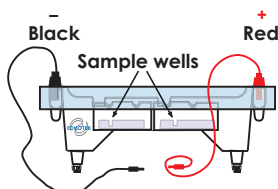

BEFORE LOADING THE SAMPLES

This experiment requires a 1.0% agarose gel and is designed for staining with InstaStain® Ethidium Bromide. Make sure the gel is completely submerged under buffer before loading the samples and conducting electrophoresis.

LOADING DNA SAMPLES

1. Load 20 μ l of the PCR sample

Lane	
1	200 bp DNA Ladder (P5)
2	Reaction sample after 30 cycles

2. Record the position of your sample in the gel for easy identification after staining.

RUNNING THE GEL

3. After the DNA samples are loaded, properly orient the cover and carefully snap it onto the electrode terminals.
4. Insert the plugs of the black and red wires into the corresponding inputs of the power source.
5. Set the power source at the required voltage and conduct electrophoresis for the length of time determined by your instructor.
6. Check to see that current is flowing properly - you should see bubbles forming on the two platinum electrodes.
7. After the electrophoresis is completed, disconnect the power and remove the gel from the bed for staining.

STAINING AND VISUALIZATION OF DNA

After electrophoresis, agarose gels require staining to visualize the separated DNA samples. Your instructor will provide instructions for DNA staining with InstaStain® Ethidium Bromide.



Duplication of this document, in conjunction with use of accompanying reagents, is permitted for classroom/laboratory use only. This document, or any part, may not be reproduced or distributed for any other purpose without the written consent of EDVOTEK, Inc. Copyright © 1996, 1997, 1998, 2000, 2001, 2002, 2005, 2007 EDVOTEK, Inc., all rights reserved

EVT 003157K

Module II: Preparation of Insert for Ligation

The amplified DNA (insert) must be cleaned (the *Taq* DNA polymerase must be removed), then cut with *Eco* RI and *Hind* III before proceeding to the ligation step.

1. Transfer the amplified DNA (~25 μ l) to a clean 1.5 ml microcentrifuge tube and label it with your group number or initials.
2. Add 300 μ l of the Purification Buffer (C1) to the tube containing the amplified DNA insert. Mix by inversion several times.
3. Place a spin column into a reservoir tube and transfer all of the mixture from step 2 to the column. Close the cap on the column.
4. Balance the tubes in a microcentrifuge and spin at maximum speed for 1 minute. At this point, the DNA product is bound to the membrane in the column.
5. Empty the reservoir tube and replace the spin column.
6. Wash the column.
 - Add 750 μ l wash buffer to the spin column.
 - Close the cap on the column.
 - Spin balanced tubes in a microcentrifuge at maximum speed for 1 minute.
7. Empty the contents of the reservoir tube and replace the spin column to the empty reservoir tube. Close the cap on the column and spin balanced tubes in a microcentrifuge at maximum speed for 1 minute.
8. Recover the DNA insert from the column:
 - Place the column into a clean 1.5 ml snap-top tube.
 - Add 12 μ l of Enzyme Grade Ultrapure Water to the center of the column (directly onto the white membrane).
 - Close the cap on the column and let the tube stand for 1 minute.
9. Balance and carefully position the tubes in the microcentrifuge so that the open caps of the snap-top tubes do not interfere with operation of the centrifuge.
10. Spin the tubes for 1 minute at maximum speed. The DNA insert will be recovered in a volume slightly less than 12 μ l.
11. Label the tube "PCR insert" and discard the used spin column.

IMPORTANT!




Make sure that the Enzyme Grade Ultrapure Water (C4) is added directly onto the white membrane of the column (avoid touching the membrane with the pipet tip) for complete elution of the DNA.



OPTIONAL STOPPING POINT

Freeze the samples at -20°C or continue with restriction enzyme digestion

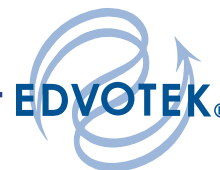

Module II: Preparation of Insert for Ligation

RESTRICTION ENZYME DIGESTION

1. To the tube labeled "PCR insert" containing 10-12 μ l the purified PCR mixture, add the following:
 - 5 μ l Restriction Enzyme Reaction Buffer
 - 10 μ l Enzyme Grade Ultrapure Water
 - 10 μ l *Eco* RI enzyme
 - 10 μ l *Hind* III enzyme
2. Mix the solution by gently pipeting up and down several times.
3. Incubate at 37°C for 60 minutes. Near the end of the incubation, prepare a separate 65°C waterbath.

After the 60 minute restriction enzyme digestion, the DNA digestion requires an incubation at 65°C to heat-inactivate the enzymes prior to being ligated into the vector.

4. Heat-inactivate the restriction enzymes by incubating the completely digested DNA at 65°C for 20 minutes.
5. Label the tube "Cut PCR Insert". After this step, the DNA insert is now ready to be ligated into the vector.



OPTIONAL STOPPING POINT

Freeze the samples at -20°C or continue with the ligation.



Duplication of this document, in conjunction with use of accompanying reagents, is permitted for classroom/laboratory use only. This document, or any part, may not be reproduced or distributed for any other purpose without the written consent of EDVOTEK, Inc. Copyright © 1996, 1997, 1998, 2000, 2001, 2002, 2005, 2007 EDVOTEK, Inc., all rights reserved
EVT 003157K

Module III: Ligation of the PCR Amplified 270 bp DNA into pUC19

1. Equilibrate an ice water bath at 16°C for Step 4.
2. In a clean microcentrifuge tube, carefully mix together:
 - 9 μ l "Cut PCR Insert" DNA
 - 3 μ l pUC19 Plasmid Vector cut with *Eco* RI and *Hind* III
 - 3 μ l Ligation buffer
 - 12 μ l Enzyme Grade Ultrapure Water
 - 3 μ l DNA Ligase30 μ l Total
Label this tube "PCR Insert+Vector/Lig".
3. In another clean microcentrifuge tube, carefully mix together:
 - 3 μ l pUC19 Plasmid Vector cut with *Eco* RI and *Hind* III
 - 3 μ l Ligation buffer
 - 21 μ l Enzyme Grade Ultrapure Water
 - 3 μ l DNA Ligase30 μ l Total
Label this tube "Vector Control/Lig".
4. Incubate the reaction tubes in a 16°C ice-water bath for 30 minutes.
Alternatively, the reaction tubes can be incubated for 2 hours at room temperature.

**OPTIONAL STOPPING POINT**

Freeze the samples at -20°C or continue with the transformation.

Module IV: Transformation

The Experiment

Note:

Sufficient colonies that have been completely suspended in CaCl_2 are critical to the success of the transformation portion of the experiment.

1. Label one microcentrifuge tube "PCR+Vector/Transf".
(This will be the transformation tube with the PCR amplified Gene.)
2. Label a second microcentrifuge tube "Vector Control/Transf".
(This will be the transformation tube with the vector control.)
3. Using a sterile pipet, add 500 μl (0.5 ml) of ice cold CaCl_2 solution to one of the tubes from step 1 or 2.
4. Pick colonies from the source plate of *E. coli* cells.
 - Use a sterile toothpick to transfer 5 colonies (2-4 mm) from the source plate to the tube containing the CaCl_2 .
 - Between your fingers, twist the tooth pick vigorously in the CaCl_2 solution to dislodge the cells.
5. Suspend the cells in the tube by tapping or vortexing (preferred).
At this point, the CaCl_2 cell suspension should look cloudy and slightly turbid. If it does not, add a few additional colonies and/or make sure the cells are completely suspended by mixing and vortexing.
6. After the cells are completely resuspended, transfer half (0.25 ml) of the cell suspension to the other labeled tube.
7. Add 10 μl "PCR+Vector/Lig" DNA to the tube labeled "PCR+Vector/Transf" and vortex.
8. Add 10 μl "Vector Control/Lig" DNA to the tube labeled "Vector Control/Transf".
9. Incubate the two tubes on ice for 20 minutes.
10. Briefly vortex the tubes, then float both tubes at 42°C for 90 seconds for the heat shock step. This facilitates the entry of DNA in bacterial cells.
11. Return both tubes immediately to ice and incubate for 2 minutes.
12. With a sterile pipet, add 250 μl of Recovery Broth to each tube & vortex.
13. Incubate cells for 30 minutes in a 37°C waterbath for a recovery period.
14. While the tubes are incubating, label 2 agar plates as follows:
 - "PCR insert + Vector"
 - "Vector Control"
 - Put your initials or group number on both plates.
15. After the recovery period, remove the tubes from the waterbath and place them on the lab bench.



Duplication of this document, in conjunction with use of accompanying reagents, is permitted for classroom/laboratory use only. This document, or any part, may not be reproduced or distributed for any other purpose without the written consent of EDVOTEK, Inc. Copyright © 1996, 1997, 1998, 2000, 2001, 2002, 2005, 2007 EDVOTEK, Inc., all rights reserved
EVT 003157K

Module IV: Transformation

Place the prepared plates inverted in a 37°C incubation oven to warm them before they are needed.

PLATING CELLS

16. Balance and centrifuge both tubes at maximum speed for 1 – 2 minutes to pellet the cells.
17. After the centrifuge has stopped, use a different pipet to carefully remove 0.4 ml of supernate from each tube. Vortex the tubes to completely resuspend the cells in the remaining liquid (0.1 ml).
18. Use a sterile pipet to transfer all of the suspended cells from the tube labeled "PCR Insert + Vector/Transf" to the middle of the plate "PCR Insert + Vector".
19. Use a sterile pipet to transfer the suspended cells from the tube labeled "Vector Control/Transf" to the middle of the plate "Vector Control".
20. Spread the cells with a sterile inoculating loop.
21. Cover both plates and allow the liquid to be absorbed.

PREPARING PLATES FOR INCUBATION

22. Stack your group's set of plates on top of one another and tape them together. Put your initials or group number on the taped set of plates.
23. Leave the plates in an upright position to allow the cell suspension to be absorbed by the agar. Place the set of plates in a safe place designated by your instructor.
24. After the cell suspension is absorbed by the agar for approximately 15-30 minutes, you or your instructor will place the plates in the inverted position (agar side on top) into a 37°C bacterial incubation oven for overnight incubation (15-20 hours).

Note: The plates are inverted to prevent condensation on the lid, which could drip onto the culture and interfere with experimental results.

AFTER OVERNIGHT INCUBATION

25. Observe the plates and estimate the number of transformants (both white and blue colonies) on each plate. Keep track of the counted colonies by putting a dot over them on the outside of the plate with a lab marker.

Module IV: Transformation

26. Calculate the transformation efficiencies for total transformants and for colonies that contain vectors with inserts (white colonies). The final volume of cells was 50 μ l and the volume plated was 50 μ l.

$$\begin{array}{rcccl} \text{Number of} & & \text{Final volume} & & \text{Number of} \\ \text{transformants} & \times & \text{of cells (ml)} & = & \text{transformants} \\ \mu\text{g of DNA} & & \text{volume plated (ml)} & & \text{per } \mu\text{g} \end{array}$$

The quantity of DNA plasmid used for the vector control was approximately 0.012 μ g. Estimate the PCR+Vector DNA to be approximately 0.05 μ g.

Study Questions

Answer the following study questions in your laboratory notebook or on a separate worksheet.

1. Why does this cloning experiment yield both blue and white colonies?
2. Do all the white and blue colonies contain a plasmid?
3. Why is there a purification step prior to the restriction enzyme digestion step?
4. Why are there two different sticky ends on the PCR amplified DNA?
5. What is a diagnostic step to analyze the success of the various steps in this experiment?



<p>EDVOTEK</p> <p>Material Safety Data Sheet</p> <p>May be used to comply with OSHA's Hazard Communication Standard, 29 CFR 1910.1200. Standard must be consulted for specific requirements.</p>		<p>IDENTITY (As Used on Label and List)</p> <p>Get loading solution concentrate, 10x</p>		<p>Back: Back exposure not permitted. If any back exposure is applicable, or no information is available, this space must be marked to indicate that.</p>	
<p>Section I</p> <p>Manufacturer's Name EDVOTEK, Inc.</p> <p>Address (Number, Street, City, State, Zip Code) 14676 Rothgeb Drive Rockville, MD 20850</p> <p>Emergency Telephone Number (301) 251-5990 Telephone Number for information (301) 251-5990 Date Prepared 10/10/06 Signature of Preparer (optional)</p>		<p>Section II - Hazardous Ingredients/Identify Information</p> <p>Hazardous Components (Specific Chemical Identity, Common Name(s)) OSHA PEL ACGIH TLV Recommended % (Optional) % (Optional)</p> <p>This product contains no hazardous materials as defined by the OSHA Hazard Communication Standard.</p>			
<p>Section III - Physical/Chemical Characteristics</p> <p>Boiling Point No data Specific Gravity (H₂O = 1) No data</p> <p>Vapor Pressure (mm Hg.) No data Melting Point N/A</p> <p>Vapor Density (AIR = 1) No data Evaporation Rate (Butyl Acetate = 1) No data</p> <p>Solubility in Water Soluble</p> <p>Appearance and Odor Blue liquid, no odor</p>		<p>Section IV - Physical/Chemical Characteristics</p> <p>Flash Point (Method Used) No data Flammable Limits No data UEL</p> <p>Extinguishing Media Dry chemical, carbon dioxide, water spray or foam</p> <p>Special Fire Fighting Procedures Use agents suitable for type of surrounding fire. Keep upwind, avoid breathing hazardous sulfur oxides and bromides. Wear SCBA.</p> <p>Unusual Fire and Explosion Hazards Unknown</p>			

<p>EDVOTEK</p> <p>Material Safety Data Sheet</p> <p>May be used to comply with OSHA's Hazard Communication Standard, 29 CFR 1910.1200. Standard must be consulted for specific requirements.</p>		<p>IDENTITY (As Used on Label and List)</p> <p>50x Electrophoresis Buffer</p>		<p>Back: Back exposure not permitted. If any back exposure is applicable, or no information is available, this space must be marked to indicate that.</p>	
<p>Section I</p> <p>Manufacturer's Name EDVOTEK, Inc.</p> <p>Address (Number, Street, City, State, Zip Code) 14676 Rothgeb Drive Rockville, MD 20850</p> <p>Emergency Telephone Number (301) 251-5990 Telephone Number for information (301) 251-5990 Date Prepared 10/05/06 Signature of Preparer (optional)</p>		<p>Section II - Hazardous Ingredients/Identify Information</p> <p>Hazardous Components (Specific Chemical Identity, Common Name(s)) OSHA PEL ACGIH TLV Recommended % (Optional) % (Optional)</p> <p>This product contains no hazardous materials as defined by the OSHA Hazard Communication Standard.</p>			
<p>Section III - Physical/Chemical Characteristics</p> <p>Boiling Point No data Specific Gravity (H₂O = 1) No data</p> <p>Vapor Pressure (mm Hg.) No data Melting Point No data</p> <p>Vapor Density (AIR = 1) No data Evaporation Rate (Butyl Acetate = 1) No data</p> <p>Solubility in Water Appreciable, (greater than 10%)</p> <p>Appearance and Odor Clear, liquid, slight vinegar odor</p>		<p>Section IV - Physical/Chemical Characteristics N.D. = No data</p> <p>Flash Point (Method Used) No data Flammable Limits LEL N.D. UEL N.D.</p> <p>Extinguishing Media Use extinguishing media appropriate for surrounding fire.</p> <p>Special Fire Fighting Procedures Wear protective equipment and SCBA with full facepiece operated in positive pressure mode.</p> <p>Unusual Fire and Explosion Hazards None identified</p>			

<p>EDVOTEK</p> <p>Material Safety Data Sheet</p> <p>May be used to comply with OSHA's Hazard Communication Standard, 29 CFR 1910.1200. Standard must be consulted for specific requirements.</p>		<p>IDENTITY (As Used on Label and List)</p> <p>Agarose</p>		<p>Back: Back exposure not permitted. If any back exposure is applicable, or no information is available, this space must be marked to indicate that.</p>	
<p>Section I</p> <p>Manufacturer's Name EDVOTEK, Inc.</p> <p>Address (Number, Street, City, State, Zip Code) 14676 Rothgeb Drive Rockville, MD 20850</p> <p>Emergency Telephone Number (301) 251-5990 Telephone Number for information (301) 251-5990 Date Prepared 10/05/06 Signature of Preparer (optional)</p>		<p>Section II - Hazardous Ingredients/Identify Information</p> <p>Hazardous Components (Specific Chemical Identity, Common Name(s)) OSHA PEL ACGIH TLV Recommended % (Optional) % (Optional)</p> <p>This product contains no hazardous materials as defined by the OSHA Hazard Communication Standard.</p>			
<p>Section III - Physical/Chemical Characteristics</p> <p>Boiling Point For 1% solution 194 F Specific Gravity (H₂O = 1) 2 No data</p> <p>Vapor Pressure (mm Hg.) No data Melting Point No data</p> <p>Vapor Density (AIR = 1) No data Evaporation Rate (Butyl Acetate = 1) No data</p> <p>Solubility in Water Insoluble - cold</p> <p>Appearance and Odor White powder, no odor</p>		<p>Section IV - Physical/Chemical Characteristics N.D. = No data</p> <p>Flash Point (Method Used) No data Flammable Limits LEL N.D. UEL N.D.</p> <p>Extinguishing Media Water spray, dry chemical, carbon dioxide, halon or standard foam</p> <p>Special Fire Fighting Procedures Possible fire hazard when exposed to heat or flame</p> <p>Unusual Fire and Explosion Hazards None</p>			

<p>Section V - Reactivity Data</p> <p>Stability Unstable Stable X Conditions to Avoid None</p> <p>Incompatibility None known</p> <p>Hazardous Decomposition or Byproducts Sulfur oxides and bromides</p> <p>Hazardous Polymerization May Occur Will Not Occur X Conditions to Avoid None</p> <p>Section VI - Health Hazard Data</p> <p>Route(s) of Entry: Inhalation? Yes Skin? Yes Ingestion? Yes</p> <p>Health Hazards (Acute and Chronic) Yes No data available for other routes</p> <p>Carcinogenicity: NTP No data IARC Monographs? OSHA Regulation? No data</p> <p>Signs and Symptoms of Exposure May cause skin or eye irritation</p> <p>Medical Conditions Generally Aggravated by Exposure None reported</p> <p>Emergency First Aid Procedures Treat symptomatically and supportively. Rinse contacted area with copious amounts of water.</p> <p>Section VII - Precautions for Safe Handling and Use</p> <p>Steps to be Taken in case Material is Released or Spilled Rinse contacted area with copious amounts of water.</p> <p>Waste Disposal Method Observe all federal, state, and local regulations.</p> <p>Precautions to be Taken in Handling and Storing Avoid eye and skin contact.</p> <p>Other Precautions None</p> <p>Section VIII - Control Measures</p> <p>Respiratory Protection (Specify Type) Chemical cartridge respirator with organic vapor cartridge.</p> <p>Ventilation Local Exhaust Yes Special Other Yes</p> <p>Protective Gloves Mechanical (General) Yes Eye Protection Other None</p> <p>Other Protective Clothing or Equipment None required</p> <p>Work/Hygiene Practices Do not ingest. Avoid contact with skin, eyes and clothing. Wash thoroughly after handling.</p>	
--	--

<p>Section V - Reactivity Data</p> <p>Stability Unstable Stable X Conditions to Avoid None</p> <p>Incompatibility Strong oxidizing agents</p> <p>Hazardous Decomposition or Byproducts Carbon monoxide, Carbon dioxide</p> <p>Hazardous Polymerization May Occur Will Not Occur X Conditions to Avoid None</p> <p>Section VI - Health Hazard Data</p> <p>Route(s) of Entry: Inhalation? Yes Skin? Yes Ingestion? Yes</p> <p>Health Hazards (Acute and Chronic) None</p> <p>Carcinogenicity: None identified NTP? IARC Monographs? OSHA Regulation?</p> <p>Signs and Symptoms of Exposure Irritation to upper respiratory tract, skin, eyes</p> <p>Medical Conditions Generally Aggravated by Exposure None</p> <p>Emergency First Aid Procedures Ingestion: If conscious, give large amounts of water. Eyes: Flush with water. Inhalation: Move to fresh air. Skin: Wash with soap and water.</p> <p>Section VII - Precautions for Safe Handling and Use</p> <p>Steps to be Taken in case Material is Released or Spilled Rinse with water, or collect in absorbent material and dispose of the absorbent material.</p> <p>Waste Disposal Method Dispose in accordance with all applicable federal, state, and local environmental regulations.</p> <p>Precautions to be Taken in Handling and Storing Avoid eye and skin contact.</p> <p>Other Precautions None</p> <p>Section VIII - Control Measures</p> <p>Respiratory Protection (Specify Type)</p> <p>Ventilation Local Exhaust Yes Special Other None</p> <p>Protective Gloves Mechanical (General) Yes Eye Protection Safety goggles</p> <p>Other Protective Clothing or Equipment None</p> <p>Work/Hygiene Practices None</p>	
--	--

<p>Section V - Reactivity Data</p> <p>Stability Unstable Stable X Conditions to Avoid None</p> <p>Incompatibility No data available</p> <p>Hazardous Decomposition or Byproducts</p> <p>Hazardous Polymerization May Occur Will Not Occur X Conditions to Avoid None</p> <p>Section VI - Health Hazard Data</p> <p>Route(s) of Entry: Inhalation? Yes Skin? Yes Ingestion? Yes</p> <p>Health Hazards (Acute and Chronic) Irritation. Large amounts may cause diarrhea</p> <p>Carcinogenicity: NTP No data IARC Monographs? OSHA Regulation?</p> <p>Signs and Symptoms of Exposure No data available</p> <p>Medical Conditions Generally Aggravated by Exposure No data available</p> <p>Emergency First Aid Procedures Treat symptomatically and supportively</p> <p>Section VII - Precautions for Safe Handling and Use</p> <p>Steps to be Taken in case Material is Released or Spilled Sweep up and place in suitable container for disposal</p> <p>Waste Disposal Method Normal solid waste disposal</p> <p>Precautions to be Taken in Handling and Storing None</p> <p>Other Precautions None</p> <p>Section VIII - Control Measures</p> <p>Respiratory Protection (Specify Type) Chemical cartridge respirator with full facepiece.</p> <p>Ventilation Local Exhaust Yes Special Other</p> <p>Protective Gloves Yes Eye Protection Splash proof goggles</p> <p>Other Protective Clothing or Equipment Impermeable clothing to prevent skin contact</p> <p>Work/Hygiene Practices None</p>	
---	--

Material Safety Data Sheets

Full size (8.5 x 11") pdf copy of MSDS available at www.edvotek.com or by request.

EDVO-Kit #
331

Material Safety Data Sheet May be used to comply with OSHA's Hazard Communication Standard, 29 CFR 1910.1200. Standard must be consulted for specific requirements.	
IDENTITY (As Used on Label and List) IPTG	Note: Blank spaces are not permitted. If any item is not applicable, or no information is available, the space must be marked to indicate that.
Section I Manufacturer's Name EDVOTEK, Inc. Address (Number, Street, City, State, Zip Code) 14676 Rothgeb Drive Rockville, MD 20850 Emergency Telephone Number (301) 251-5990 Telephone Number for information (301) 251-5990 Date Prepared 02/08/07 Signature of Preparer (optional)	Other Limits Recommended % (Optional) OSHA PEL ACGIH TLV
Section II - Hazardous Ingredients/Identify Information Hazardous Components Specific Chemical Identity, Common Name(s) OSHA PEL ACGIH TLV % (Optional) Not applicable	
Section III - Physical/Chemical Characteristics Boiling Point None Specific Gravity (H ₂ O = 1) Unknown Vapor Pressure (mm Hg) None Melting Point 109-110C Vapor Density (AIR = 1) None Evaporation Rate (Butyl Acetate = 1) None Solubility in Water Moderate	
Appearance and Odor White crystals/slight odor thioptenol	
Section IV - Physical/Chemical Characteristics Flash Point (Method Used) Unknown Flammable Limits LEL UEL Extinguishing Media Water, carbon dioxide, or dry chemical Special Fire Fighting Procedures None	
Unusual Fire and Explosion Hazards None	

Material Safety Data Sheet May be used to comply with OSHA's Hazard Communication Standard, 29 CFR 1910.1200. Standard must be consulted for specific requirements.	
IDENTITY (As Used on Label and List) Ampicillin	Note: Blank spaces are not permitted. If any item is not applicable, or no information is available, the space must be marked to indicate that.
Section I Manufacturer's Name EDVOTEK, Inc. Address (Number, Street, City, State, Zip Code) 14676 Rothgeb Drive Rockville, MD 20850 Emergency Telephone Number (301) 251-5990 Telephone Number for information (301) 251-5990 Date Prepared 02/08/07 Signature of Preparer (optional)	Other Limits Recommended % (Optional) OSHA PEL ACGIH TLV
Section II - Hazardous Ingredients/Identify Information Hazardous Components Specific Chemical Identity, Common Name(s) OSHA PEL ACGIH TLV % (Optional) Ampicillin No data	
Section III - Physical/Chemical Characteristics Boiling Point No data Specific Gravity (H ₂ O = 1) No data Vapor Pressure (mm Hg) No data Melting Point No data Vapor Density (AIR = 1) No data Evaporation Rate (Butyl Acetate = 1) No data Solubility in Water Slightly soluble	
Appearance and Odor Odorless, white crystalline powder	
Section IV - Physical/Chemical Characteristics Flash Point (Method Used) No data Flammable Limits LEL N.D. UEL N.D. Extinguishing Media Dry chemical, carbon dioxide, water spray or regular foam Special Fire Fighting Procedures Move container from fire area if possible. Do not scatter spilled material with water streams. Unusual Fire and Explosion Hazards Avoid breathing vapors.	

Material Safety Data Sheet May be used to comply with OSHA's Hazard Communication Standard, 29 CFR 1910.1200. Standard must be consulted for specific requirements.	
IDENTITY (As Used on Label and List) InstaStar® Ethidium Bromide	Note: Blank spaces are not permitted. If any item is not applicable, or no information is available, the space must be marked to indicate that.
Section I Manufacturer's Name InstaStar, Inc. P.O. Box 1232 West Bethesda, MD 20827 Emergency Telephone Number (301) 251-5990 Telephone Number for information (301) 251-5990 Date Prepared 10/05/06 Signature of Preparer (optional)	Other Limits Recommended % (Optional) OSHA PEL ACGIH TLV Data not available
Section II - Hazardous Ingredients/Identify Information Hazardous Components Specific Chemical Identity, Common Name(s) OSHA PEL ACGIH TLV % (Optional) (2,7-Diamino-10-Ethyl-9-Phenylphenanthridinium Bromide)	
Section III - Physical/Chemical Characteristics Boiling Point No data Specific Gravity (H ₂ O = 1) No data Vapor Pressure (mm Hg) No data Melting Point No data Vapor Density (AIR = 1) No data Evaporation Rate (Butyl Acetate = 1) No data Solubility in Water Soluble	
Appearance and Odor Chemical bound to paper, no odor	
Section IV - Physical/Chemical Characteristics Flash Point (Method Used) No data Flammable Limits LEL N.D. UEL N.D. Extinguishing Media Water spray, carbon dioxide, dry chemical powder, alcohol or polymer foam Special Fire Fighting Procedures Wear protective clothing and SCBA to prevent contact with skin & eyes Unusual Fire and Explosion Hazards Emits toxic fumes	

Section V - Reactivity Data Stability Unstable X Stable Conditions to Avoid None Incompatibility Strong oxidizing agents	Carbon dioxide and sulfur dioxide Conditions to Avoid None
Hazardous Decomposition or Byproducts Hazardous Polymerization May Occur X Will Not Occur Section VI - Health Hazard Data Route(s) of Entry: Inhalation? Yes Skin? Not studied Ingestion? Not studied	IARC Monographs? OSHA Regulation? No data No data
Health Hazards (Acute and Chronic) Toxicity has not been studied Carcinogenicity: Unknown Signs and Symptoms of Exposure Unknown Medical Conditions Generally Aggravated by Exposure Unknown	Emergency First Aid Procedures External: flush with water Internal: induces vomiting, consult physician
Section VII - Precautions for Safe Handling and Use Steps to be Taken in case Material is Released or Spilled Cover and sweep up with inert carrier	
Waste Disposal Method Mix material with combustible solvent and burn in a chemical incinerator with equipped afterburner and scrubber or sweep up and return to original container. Precautions to be Taken in Handling and Storing Avoid dust store cool Information CAS #367-93-1	
Section VIII - Control Measures Respiratory Protection (Specify Type) Filter mask Ventilation Local Exhaust Yes Mechanical (General) Yes Protective Gloves Rubber or vinyl Eye Protection Face mask or goggles Other Protective Clothing or Equipment Lab apron Work/Hygiene Practices Avoid dust or contact with skin	

Section V - Reactivity Data Stability Unstable X Stable Conditions to Avoid Incompatibilities Incompatibility Strong oxidizers	Toxic oxides of carbon, nitrogen and sulfur Conditions to Avoid Incompatibilities
Hazardous Decomposition or Byproducts Hazardous Polymerization May Occur X Will Not Occur Section VI - Health Hazard Data Route(s) of Entry: Inhalation? Yes Skin? Yes Ingestion? Yes	IARC Monographs? OSHA Regulation? No data No data
Health Hazards (Acute and Chronic) Sensitizers may result in allergic reaction Carcinogenicity: No data Signs and Symptoms of Exposure Repeated exposure may result in sensitization and possible anaphylactic shock. Medical Conditions Generally Aggravated by Exposure No data	Emergency First Aid Procedures Ingestion: Allergic symptoms. Eyes/Skin: Flush with water
Section VII - Precautions for Safe Handling and Use Steps to be Taken in case Material is Released or Spilled Wear suitable protective clothing. Sweep up and place in suitable container for later disposal. Do not flush spilled material down sink. Waste Disposal Method Observe all federal, state, and local regulations Precautions to be Taken in Handling and Storing Keep away from incompatible substances Other Precautions None	
Section VIII - Control Measures Respiratory Protection (Specify Type) Ventilation Local Exhaust Yes Mechanical (General) No Protective Gloves Yes Eye Protection Splash or dust proof Other Protective Clothing or Equipment Eye wash Work/Hygiene Practices Wear protective clothing and equipment to prevent contact.	

Section V - Reactivity Data Stability Unstable X Stable Conditions to Avoid None Incompatibility Strong oxidizing agents	Carbon dioxide, nitrogen oxide, hydrogen bromide gas Conditions to Avoid None
Hazardous Decomposition or Byproducts Hazardous Polymerization May Occur X Will Not Occur Section VI - Health Hazard Data Route(s) of Entry: Inhalation? Yes Skin? Yes Ingestion? Yes	IARC Monographs? OSHA Regulation? No data No data
Health Hazards (Acute and Chronic) May alter genetic material Acute: Material irritating to mucous membranes, upper respiratory tract, eyes, skin Carcinogenicity: No data available NTP?	Irritation to mucous membranes and upper respiratory tract Medical Conditions Generally Aggravated by Exposure No data
Emergency First Aid Procedures Treat symptomatically and supportively	
Section VII - Precautions for Safe Handling and Use Steps to be Taken in case Material is Released or Spilled Wear SCBA, rubber boots, rubber gloves Waste Disposal Method Mix material with combustible solvent and burn in a chemical incinerator equipped afterburner and scrubber Precautions to be Taken in Handling and Storing Use in chemical fume hood with proper protective lab gear. Other Precautions Mutagen	
Section VIII - Control Measures Respiratory Protection (Specify Type) SCBA Ventilation Local Exhaust Yes Mechanical (General) No Protective Gloves Rubber Eye Protection Chem. safety goggles Other Protective Clothing or Equipment Rubber boots Work/Hygiene Practices Use in chemical fume hood with proper protective lab gear.	

IDENTITY (As Used on Label and List)		Name, Bulk Volume and Net Weight (If you have not applicable, or no information is available, the space must be marked to indicate that)	
Enzyme Reaction Buffer		None	
Section I - Manufacturer's Information			
Manufacturer's Name	EDVOTEK, Inc.	Emergency Telephone Number	(301) 251-5990
Address (Number, Street, City, State, Zip Code)	14676 Rothgeb Drive Rockville, MD 20850	Telephone Number for information	(301) 251-5990
Date Prepared	02/08/07	Signature of Preparer (optional)	
Section II - Hazardous Ingredients/Identify Information			
Hazardous Components Specific to this Material (See Section 9)		Other Limits	
This product contains no hazardous materials as defined by the OSHA Hazard Communication Standard.		ACGIH TLV	
Section III - Physical/Chemical Characteristics			
Boiling Point	at 760mm Hg	180C	Specific Gravity (H ₂ O = 1)
Vapor Pressure (mm Hg)	20C	0.46	Melting Point
Vapor Density (AIR = 1)		2.7	Evaporation Rate (Evap/Accrate = 1)
Solubility in Water	Complete (100%)		
Appearance and Other	Colorless liquid, faint odor		
Section IV - Physical/Chemical Characteristics			
Flash Point (Method Used)	Flammable Limits	LEL	UEL
(Closed cup) 88C (192F)		3%	43%
Extinguishing Media			
Water-soluble, carbon dioxide, dry chemical, ordinary foam.			
Special Fire Fighting Procedures			
Wear SCBA with full facepiece operated in positive pressure mode. If possible, move container from fire area.			
Unusual Fire and Explosion Hazards			
Vapors may flow along surfaces to distant ignition sources and flashback. Closed containers exposed to heat may explode. Contact w/ strong oxidizers may cause fire.			

IDENTITY (As Used on Label and List)		Name, Bulk Volume and Net Weight (If you have not applicable, or no information is available, the space must be marked to indicate that)	
Enzyme Reaction Buffer		None	
Section I - Manufacturer's Information			
Manufacturer's Name	EDVOTEK, Inc.	Emergency Telephone Number	(301) 251-5990
Address (Number, Street, City, State, Zip Code)	14676 Rothgeb Drive Rockville, MD 20850	Telephone Number for information	(301) 251-5990
Date Prepared	10/10/06	Signature of Preparer (optional)	
Section II - Hazardous Ingredients/Identify Information			
Hazardous Components Specific to this Material (See Section 9)		Other Limits	
This product contains no hazardous materials as defined by the OSHA Hazard Communication Standard.		ACGIH TLV	
Section III - Physical/Chemical Characteristics			
Boiling Point	No data	Specific Gravity (H ₂ O = 1)	No data
Vapor Pressure (mm Hg)	No data	Melting Point	N/A
Vapor Density (AIR = 1)	No data	Evaporation Rate (Evap/Accrate = 1)	No data
Solubility in Water	soluble		
Appearance and Other	Clear, liquid, no odor, dry chemical, carbon dioxide, water, spray or foam.		
Section IV - Physical/Chemical Characteristics			
Flash Point (Method Used)	Flammable Limits	LEL	UEL
No data		No data	No data
Extinguishing Media			
Use extinguishing media appropriate to surrounding fire			
Special Fire Fighting Procedures			
Remove container from fire if possible.			
Unusual Fire and Explosion Hazards			
May produce toxic gases			

IDENTITY (As Used on Label and List)		Name, Bulk Volume and Net Weight (If you have not applicable, or no information is available, the space must be marked to indicate that)	
Enzyme Reaction Dilution Buffer		None	
Section I - Manufacturer's Information			
Manufacturer's Name	EDVOTEK, Inc.	Emergency Telephone Number	(301) 251-5990
Address (Number, Street, City, State, Zip Code)	14676 Rothgeb Drive Rockville, MD 20850	Telephone Number for information	(301) 251-5990
Date Prepared	11/07/06	Signature of Preparer (optional)	
Section II - Hazardous Ingredients/Identify Information			
Hazardous Components Specific to this Material (See Section 9)		Other Limits	
This product contains no hazardous materials as defined by the OSHA Hazard Communication Standard.		ACGIH TLV	
Section III - Physical/Chemical Characteristics			
Boiling Point	No data	Specific Gravity (H ₂ O = 1)	No data
Vapor Pressure (mm Hg)	No data	Melting Point	N/A
Vapor Density (AIR = 1)	No data	Evaporation Rate (Evap/Accrate = 1)	No data
Solubility in Water	soluble		
Appearance and Other	Clear, liquid, no odor, dry chemical, carbon dioxide, water, spray or foam.		
Section IV - Physical/Chemical Characteristics			
Flash Point (Method Used)	Flammable Limits	LEL	UEL
No data		No data	No data
Extinguishing Media			
Use extinguishing media appropriate to surrounding fire			
Special Fire Fighting Procedures			
Remove container from fire if possible.			
Unusual Fire and Explosion Hazards			
May produce toxic gases			

Section V - Reactivity Data		Conditions to Avoid	
Stability	Unstable	None	
	Stable	X	Heat, flame, other sources of ignition
Incompatibility			
Strong oxidizing agents, active halogen compounds, alkali metals			
Hazardous decomposition or byproduct			
Oxide of sulfur, mercuric, active halogen compounds, alkali metals			
Hazardous Polymerization			
May Occur			
Will Not Occur			
None			
Section VI - Health Hazard Data			
Route(s) of Entry:		Inhalation?	Skin? Ingestion?
		Yes	Yes
Health Hazards (Acute and Chronic)		Chronic effects:	kidney damage, liver damage
Carcinogenicity:		NTP?	OSHA Regulation?
None		No	No
Signs and Symptoms of Exposure			
Skin/eye irritation, headache, nausea, vomiting, dizziness			
Medical Conditions Generally Aggravated by Exposure			
Skin disorders			
Emergency First Aid Procedures			
Flush skin/eyes with large amounts of water. If ingested do not induce vomiting.			
Section VII - Precautions for Safe Handling and Use			
Steps to be Taken in Case Material is Released or Spilled			
Shut off ignition sources. Wear protective clothing. Use water spray to reduce vapors			
Take up with sand or other noncombustible absorbent material and dispose of properly			
Waste Disposal Method			
Dispose in accordance with all applicable federal, state, and local environmental regulations			
Precautions to be Taken in Handling and Storing			
Store in cool, dry, well ventilated flammable liquid storage area or cabinet, store above 20C			
Other Precautions			
Product may solidify at room temperature			
Section VIII - Control Measures			
Respiratory Protection (Specify Type)			
Chemical cartridge respirator w/ organic vapor cartridge			
Ventilation		Local Exhaust	Yes
		Mechanical (General)	Yes
Protective Gloves		Special	None
		Other	None
Other Protective Clothing or Equipment		Eye Protection	
		Splash Proof Safety goggles	
Work/Hygiene Practices		Avoid contact with skin/eyes	
		Keep container tightly closed.	

Section V - Reactivity Data		Conditions to Avoid	
Stability	Unstable	None	
	Stable	X	None
Incompatibility			
Copper, iron, silver salts, hydrogen peroxide, phenol, picric acid			
Formaldehyde, ether, alcohol, nitrogen oxide, strong bases, oxidizing agents.			
Toxic gases: Carbon monoxide, carbon dioxide, nitrogen oxides, chlorine, hydrogen chloride			
Hazardous Polymerization			
May Occur			
Will Not Occur			
None			
Section VI - Health Hazard Data			
Route(s) of Entry:		Inhalation?	Skin? Ingestion?
		Yes	Yes
Health Hazards (Acute and Chronic)		Toxicity has not been quantified. Sensitivity reactions (allergic) may occur by skin penetration including anaphylactic shock.	OSHA Regulation?
Carcinogenicity:		NTP?	No data
None		ARCS Monographs?	No data
Signs and Symptoms of Exposure			
Irritation to skin/eyes, mucous membranes and upper respiratory tract.			
Medical Conditions Generally Aggravated by Exposure			
Respiratory conditions			
Emergency First Aid Procedures			
Treat symptomatically and supportively			
Section VII - Precautions for Safe Handling and Use			
Steps to be Taken in Case Material is Released or Spilled			
Pick up with absorbent material. Dispose of properly.			
Waste Disposal Method			
Mix with vermiculite and dry caustic, wrap in paper and burn in a chemical incinerator equipped with afterburner and scrubber. Ignite in presence of sodium carbonate and silled lime (CaOH)			
Precautions to be Taken in Handling and Storing			
Wear protective gear to avoid skin/eye contact			
Other Precautions			
None			
Section VIII - Control Measures			
Respiratory Protection (Specify Type)			
Chemical cartridge respirator with organic vapor cartridge			
Ventilation		Local Exhaust	Yes
		Mechanical (General)	No
Protective Gloves		Special	None
		Other	None
Other Protective Clothing or Equipment		Eye Protection	
		Splash proof goggles	
Work/Hygiene Practices		Do not ingest. Avoid contact with skin, eyes and clothing. Wash thoroughly after handling.	

Section V - Reactivity Data		Conditions to Avoid	
Stability	Unstable	None	
	Stable	X	None
Incompatibility			
Copper, iron, silver salts, hydrogen peroxide, phenol, picric acid			
Formaldehyde, ether, alcohol, nitrogen oxide, strong bases, oxidizing agents.			
Toxic gases: Carbon monoxide, carbon dioxide, nitrogen oxides, chlorine, hydrogen chloride			
Hazardous Polymerization			
May Occur			
Will Not Occur			
None			
Section VI - Health Hazard Data			
Route(s) of Entry:		Inhalation?	Skin? Ingestion?
		Yes	Yes
Health Hazards (Acute and Chronic)		Toxicity has not been quantified. Sensitivity reactions (allergic) may occur by skin penetration including anaphylactic shock.	OSHA Regulation?
Carcinogenicity:		NTP?	No data
None		ARCS Monographs?	No data
Signs and Symptoms of Exposure			
Irritation to skin/eyes, mucous membranes and upper respiratory tract.			
Medical Conditions Generally Aggravated by Exposure			
Respiratory conditions			
Emergency First Aid Procedures			
Treat symptomatically and supportively			
Section VII - Precautions for Safe Handling and Use			
Steps to be Taken in Case Material is Released or Spilled			
Pick up with absorbent material. Dispose of properly.			
Waste Disposal Method			
Mix with vermiculite and dry caustic, wrap in paper and burn in a chemical incinerator equipped with afterburner and scrubber. Ignite in presence of sodium carbonate and silled lime (CaOH)			
Precautions to be Taken in Handling and Storing			
Wear protective gear to avoid skin/eye contact			
Other Precautions			
None			
Section VIII - Control Measures			
Respiratory Protection (Specify Type)			
Chemical cartridge respirator with organic vapor cartridge			
Ventilation		Local Exhaust	Yes
		Mechanical (General)	No
Protective Gloves		Special	None
		Other	None
Other Protective Clothing or Equipment		Eye Protection	
		Splash proof goggles	
Work/Hygiene Practices		Do not ingest. Avoid contact with skin, eyes and clothing. Wash thoroughly after handling.	